Oclean Smart Sonic Electric
Toothbrush User Manual

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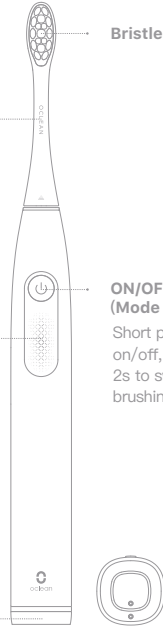
Oclean Smart Sonic Electric Toothbrush User Manual

Brush head
When replacing brush head, pull out or insert it vertically and do not rotate.

Mode Indicator
Level 1—Sensitive Cleaning
Level 2—Whitening
Level 3—Strong Cleaning

Magnetic charging base
2 hours for full charging and stand-by about 30 days.

Instructions

- 1. Power On**
Start the toothbrush by short pressing the power button, the toothbrush will automatically shut down after a 2-mins brushing plan is finished.
- 2. Change Brushing Area**
When the toothbrush is working, there will be a reminder with vibration to remind the user to change the brushing area every 30 seconds.
- 3. Intensity Switch**
Oclean F1 has preset 3 types of 2mins brushing plan: Sensitive Cleaning(Default), Whitening Cleaning, Strong Cleaning.
To switch the intensity, just long press the power button for 2s under working condition.
- 4. Battery Indication**
If battery power is less than 5%, the red indicator light will flash as a charging reminder. The indicator flashes when battery percentage is lower than 20%, needs to be charged to ensure normal use.

Low battery reminder: The toothbrush will not work when battery is lower than 5%, short press on the power button, red light flashes twice and motor vibrates twice to remind that battery is extremely low, recharge is needed in this case.

The red indicator flashes when battery is low. The white indicator is on constantly when battery is full.

Specifications

Product Name	Oclean Sonic Electric Toothbrush
Model	Oclean F1
Max Torsion	140 gf.cm
Waterproof	IPX7
Power	1.8W
Charger	USB charger
Input	5V \square 0.6A
Battery Life*	30 days
Battery Capacity	800 mAh
Charging Time	2 hours

* Calculated by 2 minutes/time, twice/day.

Facebook Page

www.facebook.com/ocleanofficial

Follow the Oclean Official Facebook Page

Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
- × The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
- * Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

Name of parts	Toxic and harmful substance or element						Environment Protection Term
	Pb	Hg	Cd	Cr6+	PBB	PBDE	
Brush head	○	○	○	○	○	○	Ten Years
Handle	×	○	○	○	○	○	
Charging base	×	○	○	○	○	○	Five Years
Accessory	×	○	○	○	○	○	
Battery	○	○	○	○	○	○	

Warranty Card

- These clauses only apply to Oclean products sold in the country or region where the purchase is made, if there is any inconsistency or omission with the national policy, the national policy shall prevail.
- For performance issues caused by the product itself, one may avail service as per below regulations with valid proof of purchase:
- If a customer can avail warranty service for free within 24 months from the date of purchase if there is a performance failure caused by the product itself.
- In order to avail replacement and warranty services, one should show the valid proof of purchase.
- Other product related materials (such as Packaging box, Product Manual etc.) are not covered in warranty.

Product model:	Serial Number:					
	Name:					
User Information	Tel:					
	Email:					
	Address:					
	Date of purchase:					
Problem Description:	Invoice number:					

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Specifications

Product Name	Oclean Sonic Electric Toothbrush
Model	Oclean F1
Max Torsion	140 gf.cm
Waterproof	IPX7
Power	1.8W
Charger	USB charger
Input	5V = 0.6A
Battery Life*	30 days
Battery Capacity	800 mAh
Charging Time	2 hours

* Calculated by 2 minutes/time, twice/day.

Facebook Page

www.facebook.com/ocleanofficial

Follow the Oclean Official Facebook Page

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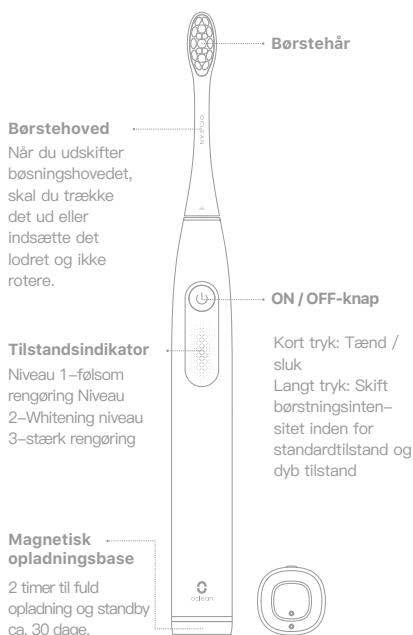
Specifications

Product Name	Oclean Sonic Electric Toothbrush
Model	Oclean F1
Max Torsion	140 gf.cm
Waterproof	IPX7
Power	1.8W
Charger	USB charger
Input	5V \square 0.6A
Battery Life*	30 days
Battery Capacity	800 mAh
Charging Time	2 hours

* Calculated by 2 minutes/time, twice/day.

Facebook Page

www.facebook.com/ocleanofficial

Follow the Oclean Official Facebook Page

Name of parts	Toxic and harmful substance or element						Environment Protection Term
	Pb	Hg	Cd	Cr6+	PBB	PBDE	
Brush head	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	Ten Years
Handle	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	
Charging base	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	
Accessory	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	
Battery	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	

Remark

- The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;
- × The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;
- * Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.

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Specifications

Product Name	Oclean Sonic Electric Toothbrush
Model	Oclean F1
Max Torsion	140 gf.cm
Waterproof	IPX7
Power	1.8W
Charger	USB charger
Input	5V \equiv 0.6A
Battery Life*	30 days
Battery Capacity	800 mAh
Charging Time	2 hours

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Facebook Page

www.facebook.com/ocleanofficial

Follow the OcClean Official Facebook Page

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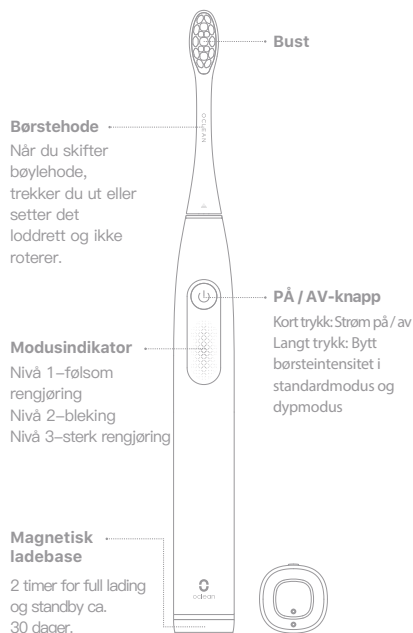
Specifications

Product Name	Oclean Sonic Electric Toothbrush
Model	Oclean F1
Max Torision	140 gf.cm
Waterproof	IPX7
Power	1.8W
Charger	USB charger
Input	5V --- 0.6A
Battery Life*	30 days
Battery Capacity	800 mAh
Charging Time	2 hours

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